


Anti-TRAF4 Rabbit pAb




GB111015 -20°C


Product Information			
Protein full name	TNF receptor-associated factor 4		
Synonyms	Cysteine-rich motif associated to RING and Traf domains protein 1, Traf4, Cart1, MLN62, tumor necrosis receptor-associated factor 4A		
Immunogen	KLH conjugated Synthetic peptide corresponding to Human TRAF4		
Uniprot ID	Q9BUZ4		
Isotype	IgG		
Purity	Affinity purification		
Subcellular location	Nucleus, Cytoplasm		
Predicted MW./Observed MW.	53 kDa / 53 kDa		
Applications	Species	Dilution	Positive Tissue
WB	Human, Mouse, Rat	1: 500-1: 1000	thymus, brain, lymph node
IHC/IF	Mouse, Rat	1: 150-1: 500	thymus, lymphatic, lung
Storage			
Storage	Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		
Storage Buffer	PBS with 0.15% ProClin300, 100 µg/mL BSA and 50% glycerol.		
background			
The tumor necrosis factor receptor-associated factor (TRAF) family is a group of signaling adaptors which participate in various signaling pathways through binding of the tumor necrosis factor receptor super family. As one of the important member of TRAF family proteins, TRAF4 was initially isolated from breast carcinomas and identified as the first member of the TRAF family to be up-regulated in human carcinomas.			
Images			

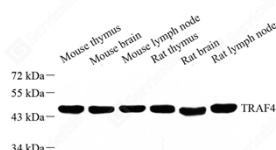

WB analysis of TRAF4(GB111015).
Sample: Protein treated by RIPA Lysis Buffer (G2002).
Blocking buffer: 3% Nonfat dry milk (GC310001) in TBST, RT, 1h.
Primary antibody: 1: 1000, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 5000, RT, 1h.

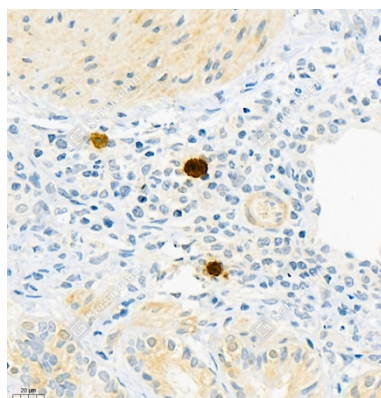

IHC analysis of TRAF4(GB111015).
Sample: Rat lung, 4% PFA (G1101) 12-24h.
Antigen retrieval: Citrate buffer (pH 6.0) (G1201), 98°C, 20 min.
Primary antibody: 1: 250, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 200, RT, 1h.

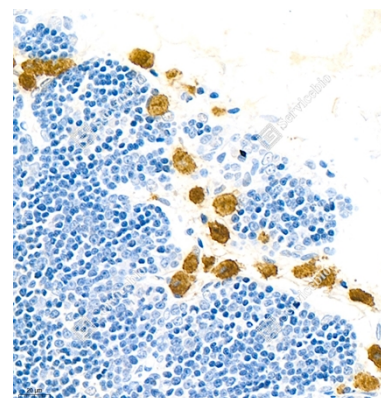

IHC analysis of TRAF4(GB111015).
Sample: Rat thymus, 4% PFA (G1101) 12-24h.
Antigen retrieval: Citrate buffer (pH 6.0) (G1201), 98°C, 20 min.
Primary antibody: 1: 250, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 200, RT, 1h.

For Research Use Only!

Ver. No.: V2.0-20260420