


Anti-BMPR2 Rabbit pAb



GB111333 -20°C

Product Information			
Protein full name	Bone morphogenetic protein receptor type-2		
Synonyms	BMP type-2 receptor, BMPR-2, BRK-3, Bone morphogenetic protein receptor type II, BMP type II receptor, BMPR2, BMPR-II		
Immunogen	Recombinant protein corresponding to Mouse BMPR2		
Uniprot ID	O35607		
Isotype	IgG		
Purity	Affinity purification		
Predicted MW./Observed MW.	115 kDa / 150 kDa		
Applications	Species	Dilution	Positive Tissue
WB	Mouse, Rat	1: 500-1: 1000	lung, brain
Storage			
Storage	Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		
Storage Buffer	PBS with 0.15% ProClin300, 100 µg/mL BSA and 50% glycerol.		
background			
BMPs (bone morphogenetic protein) are involved in endochondral bone formation and embryogenesis. BMPR2 encodes a member of the (BMP) receptor family of transmembrane serine/threonine kinases. It is a widely expressed receptor, with high mRNA expression during development in many tissues.			
Images			

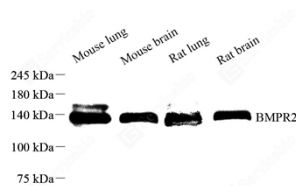

WB analysis of BMPR2(GB111333).
Sample: Protein treated by RIPA Lysis Buffer (G2002).
Blocking buffer: 3% Nonfat dry milk (GC310001) in TBST, RT, 1h.
Primary antibody: 1: 500, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 5000, RT, 1h.

For Research Use Only!

Ver. No.: V2.0-20260420