


Anti-AZI2 Rabbit pAb




GB111682 -20°C


Product Information			
Protein full name	5-azacytidine-induced protein 2		
Synonyms	NF-kappa-B-activating kinase-associated protein 1, Nak-associated protein 1, Azi2, Az2, Nap1, Tbkp2, AA410145 , MGC112644		
Immunogen	Recombinant protein corresponding to Mouse AZI2		
Uniprot ID	Q9QYP6		
Isotype	IgG		
Purity	Affinity purification		
Subcellular location	Cytoplasm		
Applications	Species	Dilution	Positive Tissue
IHC/IF	Mouse	1: 400-1: 1300	brain, kidney, testis
Storage			
Storage	Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		
Storage Buffer	PBS with 0.15% ProClin300, 100 µg/mL BSA and 50% glycerol.		
background			
AZI2, also known as NAP1, is a TNF receptor (TNFR)-associated factor family member-associated NF-κB activator-binding kinase 1-binding protein that regulates the production of IFNs. It is an adapter protein which binds TBK1 and IKBKE playing a role in antiviral innate immunity.			
Images			

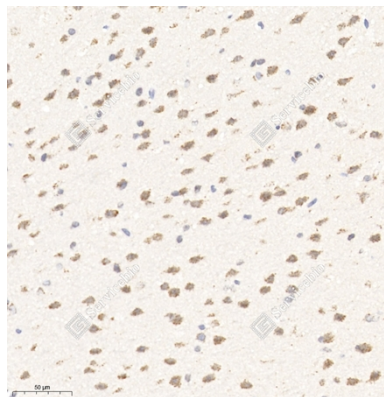

IHC analysis of AZI2(GB111682).
Sample: Mouse brain, 4% PFA (G1101) 12-24h.
Antigen retrieval: Citrate buffer (pH 6.0) (G1201), 98°C, 20 min.
Primary antibody: 1: 650, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 200, RT, 1h.

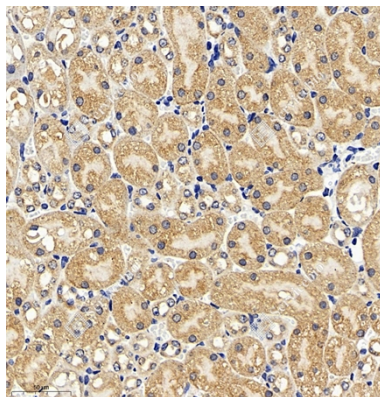

IHC analysis of AZI2(GB111682).
Sample: Mouse kidney, 4% PFA (G1101) 12-24h.
Antigen retrieval: Citrate buffer (pH 6.0) (G1201), 98°C, 20 min.
Primary antibody: 1: 650, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 200, RT, 1h.

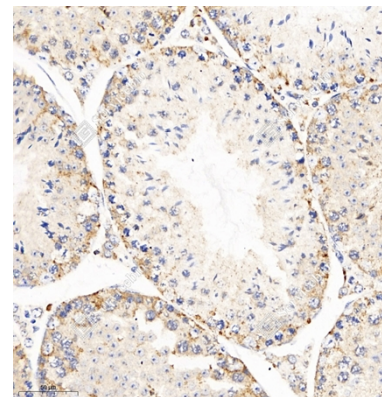

IHC analysis of AZI2(GB111682).
Sample: Mouse testis, 4% PFA (G1101) 12-24h.
Antigen retrieval: Citrate buffer (pH 6.0) (G1201), 98°C, 20 min.
Primary antibody: 1: 650, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 200, RT, 1h.

For Research Use Only!

Ver. No.: V2.0-20260420