

Anti-TRAF6 Rabbit pAb




GB112201 -20°C


Product Information			
Protein full name	TNF receptor-associated factor 6		
Synonyms	E3 ubiquitin-protein ligase TRAF6, RING-type E3 ubiquitin transferase TRAF6, Traf6, RNF85, RING finger protein 85, Interleukin-1 signal transducer		
Immunogen	Recombinant protein corresponding to Mouse TRAF6		
Uniprot ID	P70196		
Isotype	IgG		
Purity	Affinity purification		
Subcellular location	Cytoplasm, Nucleus		
Predicted MW./Observed MW.	60 kDa / 60 kDa		
Applications	Species	Dilution	Positive Tissue
WB	Human, Mouse, Rat	1: 500-1: 1000	liver, brain
IHC/IF	Mouse, Rat	1: 400-1: 800	lung, substantia nigra
Storage			
Storage	Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		
Storage Buffer	PBS with 0.15% ProClin300, 100 µg/mL BSA and 50% glycerol.		
background			
The protein encoded by this gene is a member of the TNF receptor associated factor (TRAF) protein family. TRAF proteins are associated with, and mediate signal transduction from members of the TNF receptor superfamily. This protein mediates the signaling not only from the members of the TNF receptor superfamily, but also from the members of the Toll/IL-1 family.			
Images			



WB analysis of TRAF6(GB112201).
Sample: Protein treated by RIPA Lysis Buffer (G2002).
Blocking buffer: 3% Nonfat dry milk (GC310001) in TBST, RT, 1h.
Primary antibody: 1: 500, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 5000, RT, 1h.

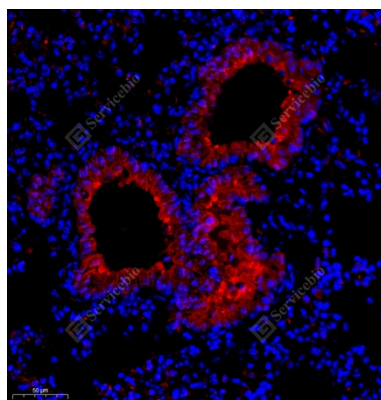

IF analysis of TRAF6(GB112201)red.
Sample: Mouse lung, 4% PFA (G1101) 12-24h.
Antigen retrieval: Tris-EDTA buffer (pH 9.0) (G1203), 98°C, 20min.
Blocking buffer: 3% BSA (GC305010) in PBS, RT, 30min.
Primary antibody: 1: 400, 4°C overnight.
Secondary antibody: Cy3 conjugated Goat Anti-Rabbit IgG (H+L) (GB21303), 1: 300, RT, 1h.

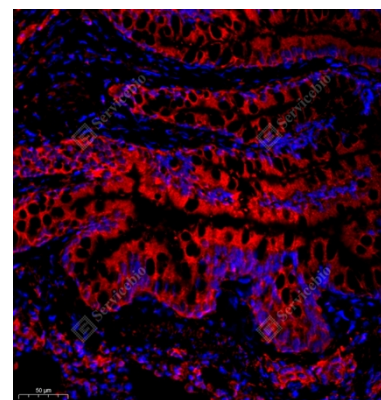

IF analysis of TRAF6(GB112201)red.
Sample: Rat lung, 4% PFA (G1101) 12-24h.
Antigen retrieval: Tris-EDTA buffer (pH 9.0) (G1203), 98°C, 20min.
Blocking buffer: 3% BSA (GC305010) in PBS, RT, 30min.
Primary antibody: 1: 400, 4°C overnight.
Secondary antibody: Cy3 conjugated Goat Anti-Rabbit IgG (H+L) (GB21303), 1: 300, RT, 1h.

For Research Use Only!

Ver. No.: V2.0-20260420