


Anti-TRAF3IP3 Rabbit pAb



GB112228 -20°C

Product Information			
Protein full name	TRAF3-interacting JNK-activating modulator		
Synonyms	TRAF3-interacting protein 3, Traf3ip3, T3jam		
Immunogen	Recombinant protein corresponding to Mouse TRAF3IP3		
Uniprot ID	Q8C0G2		
Isotype	IgG		
Purity	Affinity purification		
Predicted MW./Observed MW.	59 kDa / 59 kDa		
Applications	Species	Dilution	Positive Tissue
WB	Human, Mouse, Rat	1: 600-1: 1000	spleen, thymus, lymph node
Storage			
Storage	Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		
Storage Buffer	PBS with 0.15% ProClin300, 100 µg/mL BSA and 50% glycerol.		
background			
The gene encodes a protein that mediates cell growth by modulating the c-Jun N-terminal kinase signal transduction pathway. The encoded protein may also interact with a large multi-protein assembly containing the phosphatase 2A catalytic subunit. Alternative splicing results in multiple transcript variants.			
Images			

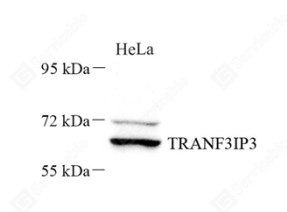

WB analysis of TRAF3IP3(GB112228).
Sample: Protein treated by RIPA Lysis Buffer (G2002).
Blocking buffer: 3% Nonfat dry milk (GC310001) in TBST, RT, 1h.
Primary antibody: 1: 600, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 5000, RT, 1h.



WB analysis of TRAF3IP3(GB112228).
Sample: Protein treated by RIPA Lysis Buffer (G2002).
Blocking buffer: 3% Nonfat dry milk (GC310001) in TBST, RT, 1h.
Primary antibody: 1: 600, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 5000, RT, 1h.

For Research Use Only!

Ver. No.: V2.0-20260420