


Anti-BinCARD Rabbit pAb



GB112417 -20°C

Product Information			
Protein full name	Caspase recruitment domain-containing protein 19		
Synonyms	Bcl10-interacting CARD protein, BinCARD, Card19, C9orf89, Chromosome 9 open reading frame 89		
Immunogen	Recombinant protein corresponding to Mouse BinCARD		
Uniprot ID	Q9D1I2		
Isotype	IgG		
Purity	Affinity purification		
Predicted MW./Observed MW.	21 kDa / 21 kDa		
Applications	Species	Dilution	Positive Tissue
WB	Mouse, Rat	1: 300-1: 600	kidney
Storage			
Storage	Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		
Storage Buffer	PBS with 0.15% ProClin300, 100 µg/mL BSA and 50% glycerol.		
background			
This gene plays a role in inhibiting the effects of BCL10-induced activation of NF-kappa-B.			
Images			

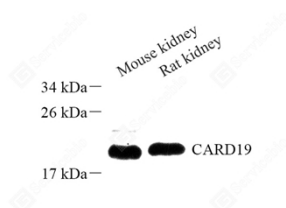

WB analysis of BinCARD(GB112417).
 Sample: Protein treated by RIPA Lysis Buffer (G2002).
 Blocking buffer: 3% Nonfat dry milk (GC310001) in TBST, RT, 1h.
 Primary antibody: 1: 300, 4°C overnight.
 Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 5000, RT, 1h.

For Research Use Only!

Ver. No.: V2.0-20260420