



Anti-URP2/Kindlin-3 Rabbit pAb

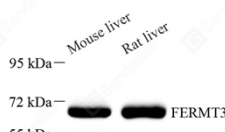

GB113250 -20°C

Product Information			
Protein full name	Fermitin family homolog 3		
Synonyms	FERMT3, KIND3, Kindlin 3, MIG 2, MIG2 like protein, MIG2B, Unc 112 related protein 2, UNC112C, URP2, URP2SF		
Immunogen	Recombinant protein corresponding to Mouse URP2/Kindlin-3		
Uniprot ID	Q8K1B8		
Isotype	IgG		
Purity	Affinity purification		
Predicted MW./Observed MW.	76 kDa / 72 kDa		
Applications	Species	Dilution	Positive Tissue
WB	Human, Mouse, Rat	1: 500-1: 1000	liver
Storage			
Storage	Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		
Storage Buffer	PBS with 0.15% ProClin300, 100 µg/mL BSA and 50% glycerol.		
background			
Kindlin-3 (FERMT-3) is known to be central in hemostasis and thrombosis control and its deficiency disrupts platelet aggregation and causes Leukocyte Adhesion Deficiency disease. Point mutations in the kindlin-3 gene have been identified in humans with a rare inherited Integrin Activation efficiency Disease (IADD), also designated as Leukocyte Adhesion Deficiency syndrome (LAD-III or LAD1 variant). Kindlin-3 expression is restricted to cells of hematopoietic origin.			
Images			

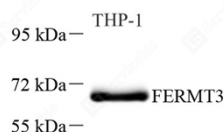

WB analysis of URP2/Kindlin-3(GB113250).
Sample: Protein treated by RIPA Lysis Buffer (G2002).
Blocking buffer: 3% Nonfat dry milk (GC310001) in TBST, RT, 1h.
Primary antibody: 1: 1000, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 5000, RT, 1h.

WB analysis of URP2/Kindlin-3(GB113250).
Sample: Protein treated by RIPA Lysis Buffer (G2002).
Blocking buffer: 3% Nonfat dry milk (GC310001) in TBST, RT, 1h.
Primary antibody: 1: 1000, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 5000, RT, 1h.

For Research Use Only!

Ver. No.: V2.0-20260420