


Anti-Rel B Rabbit pAb




GB11353 -20°C


Product Information			
Protein full name	Transcription factor RelB		
Synonyms	Reticuloendotheliosis viral oncogene homolog B, Nuclear factor of kappa light polypeptide gene enhancer in B cells 3, IREL, v rel avian reticuloendotheliosis viral oncogene homolog, RELB, I-REL, IREL, REL-B, RELB proto-oncogene, NF-kB subunit		
Immunogen	KLH conjugated Synthetic peptide corresponding to Mouse Rel B		
Uniprot ID	Q04863		
Isotype	IgG		
Purity	Affinity purification		
Subcellular location	Nucleus, Cytoplasm		
Applications	Species	Dilution	Positive Tissue
IHC/IF	Mouse, Rat	1: 1000-1: 2000	brain, colon, liver
Storage			
Storage	Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		
Storage Buffer	PBS with 0.15% ProClin300, 100 µg/mL BSA and 50% glycerol.		
background			
RelB is a pleiotropic transcription factor which is present in almost all cell types and is involved in many biological processes such as inflammation, immunity, differentiation, cell growth, tumorigenesis and apoptosis. NF- kappa-B is a homo- or heterodimeric complex formed by the Rel-like domain-containing proteins RELA/p65, RELB, NFKB1/p105, NFKB1/p50, REL and NFKB2/p52. The dimers bind at kappa-B sites in the DNA of their target genes and the individual dimers have distinct preferences for different kappa-B sites that they can bind with distinguishable affinity and specificity.			
Images			

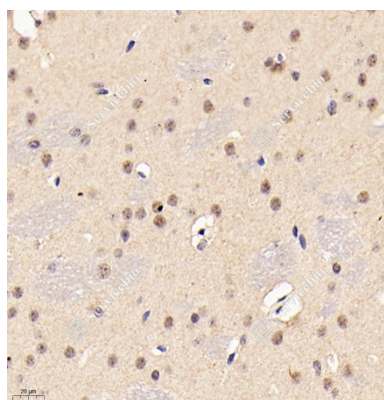

IHC analysis of Rel B(GB11353).
Sample: Mouse brain, 4% PFA (G1101) 12-24h.
Antigen retrieval: Citrate buffer (pH 6.0) (G1201), 98°C, 20 min.
Primary antibody: 1: 1000, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 200, RT, 1h.

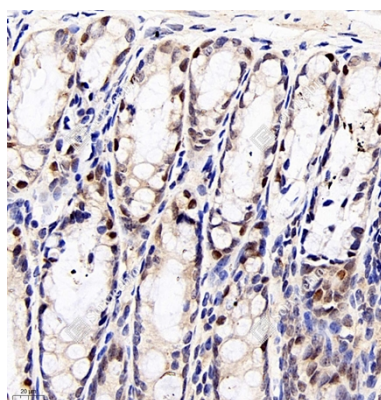

IHC analysis of Rel B(GB11353).
Sample: Mouse colon, 4% PFA (G1101) 12-24h.
Antigen retrieval: Citrate buffer (pH 6.0) (G1201), 98°C, 20 min.
Primary antibody: 1: 1000, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 200, RT, 1h.

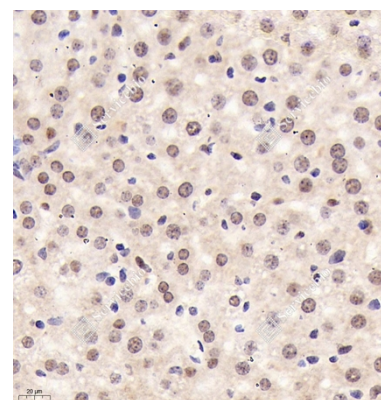

IHC analysis of Rel B(GB11353).
Sample: Mouse liver, 4% PFA (G1101) 12-24h.
Antigen retrieval: Citrate buffer (pH 6.0) (G1201), 98°C, 20 min.
Primary antibody: 1: 1000, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 200, RT, 1h.

For Research Use Only!

Ver. No.: V2.0-20260420