


Anti-TRK fused gene Rabbit pAb



GB113643 -20°C

Product Information			
Protein full name	Protein TFG		
Synonyms	FLJ36137, TF6, TFG, TRK fused gene, TRK-fused gene protein, TRKT3		
Immunogen	Recombinant protein corresponding to Human TRK fused gene		
Uniprot ID	Q92734		
Isotype	IgG		
Purity	Affinity purification		
Predicted MW./Observed MW.	43 kDa / 55 kDa		
Applications	Species	Dilution	Positive Tissue
WB	Human	1: 500-1: 1000	293T, HepG2
Storage			
Storage	Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		
Storage Buffer	PBS with 0.15% ProClin300, 100 µg/mL BSA and 50% glycerol.		
background			
Protein TFG (TRK-fused gene protein) plays a role in regulating phosphotyrosine-specific phosphatase-1 activity. Mutations in TFG may have important clinical relevance for current therapeutic strategies to treat metastatic melanoma. Defects in TFG are a cause of thyroid papillary carcinoma (TPC), a common tumor of the thyroid that typically arises as an irregular, solid or cystic mass from otherwise normal thyroid tissue.			
Images			

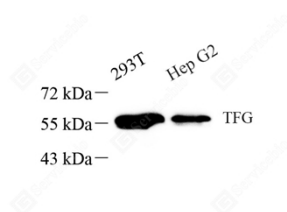

WB analysis of TRK fused gene(GB113643).
Sample: Protein treated by RIPA Lysis Buffer (G2002).
Blocking buffer: 3% Nonfat dry milk (GC310001) in TBST, RT, 1h.
Primary antibody: 1: 1000, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 5000, RT, 1h.

For Research Use Only!

Ver. No.: V2.0-20260420