


Anti-TRAF7 Rabbit pAb



GB114510 -20°C

Product Information			
Protein full name	E3 ubiquitin-protein ligase TRAF7		
Synonyms	RFWD1, RING finger protein 119, RNF119, TRAF7		
Immunogen	Recombinant protein corresponding to Human TRAF7		
Uniprot ID	Q6Q0C0		
Isotype	IgG		
Purity	Affinity purification		
Predicted MW./Observed MW.	75 kDa / 67 kDa		
Applications	Species	Dilution	Positive Tissue
WB	Mouse	1: 300-1: 600	colon
Storage			
Storage	Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		
Storage Buffer	PBS with 0.15% ProClin300, 100 µg/mL BSA and 50% glycerol.		
background			
Plays a role in the phosphorylation of MAPK1 and/or MAPK3, probably via its interaction with MAP3K3.			
Images			

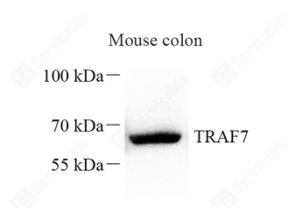

WB analysis of TRAF7(GB114510).
 Sample: Protein treated by RIPA Lysis Buffer (G2002).
 Blocking buffer: 3% Nonfat dry milk (GC310001) in TBST, RT, 1h.
 Primary antibody: 1: 600, 4°C overnight.
 Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 5000, RT, 1h.

For Research Use Only!

Ver. No.: V2.0-20260420