


Anti-NFKBID Rabbit pAb



GB114594 -20°C

Product Information			
Protein full name	NF-kappa-B inhibitor delta		
Synonyms	I-kappa-B-delta, IkappaBdelta, IkappaBNS, Ikb-delta, IKBNS, NF-kappa-B inhibitor delta, NFKBID, TA-NFKBH, T-cell activation NFKB-like protein		
Immunogen	KLH conjugated Synthetic peptide corresponding to Mouse NFKBID		
Uniprot ID	Q2TB02		
Isotype	IgG		
Purity	Affinity purification		
Predicted MW./Observed MW.	35 kDa / 35 kDa		
Applications	Species	Dilution	Positive Tissue
WB	Mouse, Rat	1: 300-1: 800	lung
Storage			
Storage	Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		
Storage Buffer	PBS with 0.15% ProClin300, 100 µg/mL BSA and 50% glycerol.		
background			
Regulates the expression of IL-2, IL-6, and other cytokines through regulation on NF-kappa-B activity. Functions in the regulation of inflammatory responses . Involved in the induction of T helper 17 cells (Th17) differentiation upon recognition of antigen by T cell antigen receptor (TCR) .			
Images			

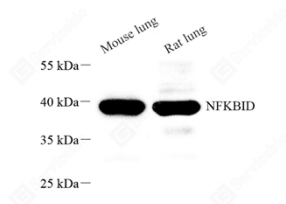

WB analysis of NFKBID(GB114594).
Sample: Protein treated by RIPA Lysis Buffer (G2002).
Blocking buffer: 3% Nonfat dry milk (GC310001) in TBST, RT, 1h.
Primary antibody: 1: 800, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 5000, RT, 1h.

For Research Use Only!

Ver. No.: V2.0-20260420