

Anti-NKRF Rabbit pAb




GB114665 -20°C

Product Information			
Protein full name	NF-kappa-B-repressing factor		
Synonyms	ITBA4, NF kappa B repressing factor, NKRF, NRF, Protein ITBA4, Transcription factor NRF		
Immunogen	KLH conjugated Synthetic peptide corresponding to Mouse NKRF		
Uniprot ID	Q8BY02		
Isotype	IgG		
Purity	Affinity purification		
Predicted MW./Observed MW.	78 kDa / 100 kDa		
Applications	Species	Dilution	Positive Tissue
WB	Human, Mouse, Rat	1: 300-1: 1000	brain
Storage			
Storage	Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		
Storage Buffer	PBS with 0.15% ProClin300, 100 µg/mL BSA and 50% glycerol.		
background			
Interacts with a specific negative regulatory element (NRE) 5'-AATTCCTCTGA-3' to mediate transcriptional repression of certain NK-kappa-B responsive genes. Involved in the constitutive silencing of the interferon beta promoter, independently of the virus-induced signals, and in the inhibition of the basal and cytokine-induced iNOS promoter activity. Also involved in the regulation of IL-8 transcription .			
Images			



WB analysis of NKRF(GB114665).
Sample: Protein treated by RIPA Lysis Buffer (G2002).
Blocking buffer: 3% Nonfat dry milk (GC310001) in TBST, RT, 1h.
Primary antibody: 1: 1000, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 5000, RT, 1h.

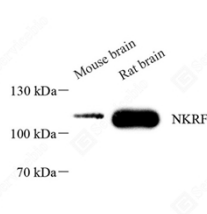

WB analysis of NKRF(GB114665).
Sample: Protein treated by RIPA Lysis Buffer (G2002).
Blocking buffer: 3% Nonfat dry milk (GC310001) in TBST, RT, 1h.
Primary antibody: 1: 1000, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 5000, RT, 1h.

For Research Use Only!

Ver. No.: V2.0-20260420