


Anti-Aggrecan Rabbit pAb




GB115687 -20°C


Product Information			
Protein full name	Aggrecan core protein		
Synonyms	ACAN, AGC1, AGCAN, CSPG1, CSPGCP, MSK16, SEDK, Aggrecan		
Immunogen	Recombinant protein corresponding to Mouse Aggrecan		
Uniprot ID	Q61282		
Isotype	IgG		
Purity	Affinity purification		
Subcellular location	Secreted, extracellular space, extracellular matrix		
Applications	Species	Dilution	Positive Tissue
IHC/IF	Human, Mouse, Rat	1: 500-1: 1000	cartilage, bone
Storage			
Storage	Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		
Storage Buffer	PBS with 0.15% ProClin300, 100 µg/mL BSA and 50% glycerol.		
background			
Aggrecan (ACAN), also known as cartilage-specific proteoglycan core protein (CSPCP) or chondroitin sulfate proteoglycan 1, is a protein that in humans is encoded by the ACAN gene. This gene is a member of the lectican (chondroitin sulfate proteoglycan) family. This proteoglycan is a major component of extracellular matrix of cartilagenous tissues. A major function of this protein is to resist compression in cartilage. It binds avidly to hyaluronic acid via an N-terminal globular region. May play a regulatory role in the matrix assembly of the cartilage.			
Images			

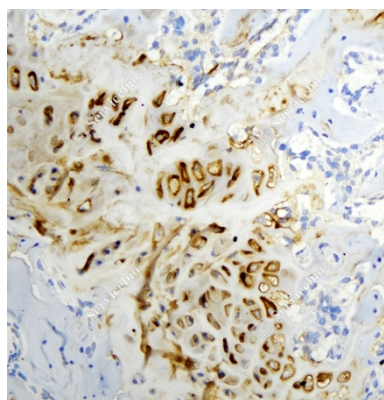

IHC analysis of Aggrecan(GB115687).
Sample: Mouse cartilage, 4% PFA (G1101) 12-24h.
Antigen retrieval: Citrate buffer (pH 6.0) (G1201), 98°C, 20 min.
Primary antibody: 1: 500, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 200, RT, 1h.

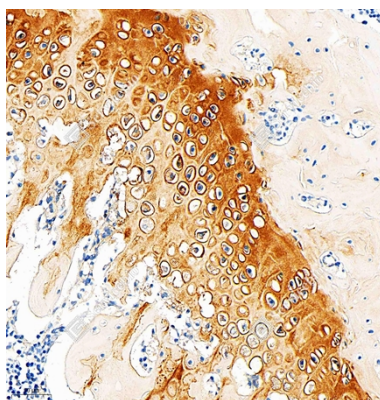

IHC analysis of Aggrecan(GB115687).
Sample: Rat cartilage, 4% PFA (G1101) 12-24h.
Antigen retrieval: Citrate buffer (pH 6.0) (G1201), 98°C, 20 min.
Primary antibody: 1: 1000, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 200, RT, 1h.

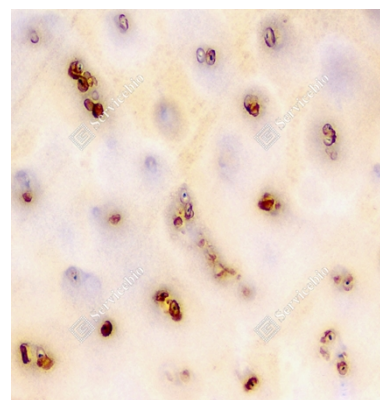

IHC analysis of Aggrecan(GB115687).
Sample: Human bone, 4% PFA (G1101) 12-24h.
Antigen retrieval: Citrate buffer (pH 6.0) (G1201), 98°C, 20 min.
Primary antibody: 1: 1000, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 200, RT, 1h.

For Research Use Only!

Ver. No.: V2.0-20260420