

Anti-RANKL Rabbit pAb



GB115705 -20°C

Product Information

Protein full name	Tumor necrosis factor ligand superfamily member 11
Synonyms	CD254, hRANKL2, ODF, OPGL, OPTB2, Osteoprotegerin ligand, RANKL, sOdf, TNFSF11, TRANCE, Opgl, Rankl, Trance
Immunogen	Recombinant protein corresponding to Mouse RANKL
Uniprot ID	O35235
Isotype	IgG
Purity	Affinity purification
Subcellular location	Cell membrane, Cytoplasm, Secreted

Applications	Species	Dilution	Positive Tissue
IHC/IF	Mouse, Rat	1: 600-1: 2400	bone, spleen


Storage

Storage	Store at -20°C for one year. Avoid repeated freeze/thaw cycles.
Storage Buffer	PBS with 0.15% ProClin300, 100 µg/mL BSA and 50% glycerol.


background

Cytokine that binds to TNFRSF11B/OPG and to TNFRSF11A/RANK. Osteoclast differentiation and activation factor. Augments the ability of dendritic cells to stimulate naive T-cell proliferation. May be an important regulator of interactions between T-cells and dendritic cells and may play a role in the regulation of the T-cell-dependent immune response. May also play an important role in enhanced bone-resorption in humoral hypercalcemia of malignancy. Induces osteoclastogenesis by activating multiple signaling pathways in osteoclast precursor cells, chief among which is induction of long lasting oscillations in the intracellular concentration of Ca²⁺ resulting in the activation of NFATC1, which translocates to the nucleus and induces osteoclast-specific gene transcription to allow differentiation of osteoclasts. During osteoclast differentiation, in a TMEM64 and ATP2A2-dependent manner induces activation of CREB1 and mitochondrial ROS generation necessary for proper osteoclast generation.


Images

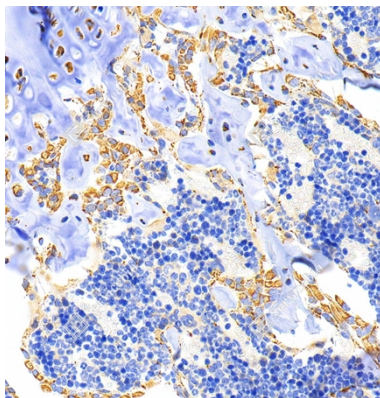

IHC analysis of RANKL(GB115705).
Sample: Mouse bone, 4% PFA (G1101) 12-24h.
Antigen retrieval: Citrate buffer (pH 6.0) (G1201), pressure cooker 2 min.
Primary antibody: 1: 1000, 4°C overnight.
Secondary antibody: S-vision poly-HRP conjugated Goat Anti-Rabbit IgG (H+L), Ready to use(G1302), RT, 20min.

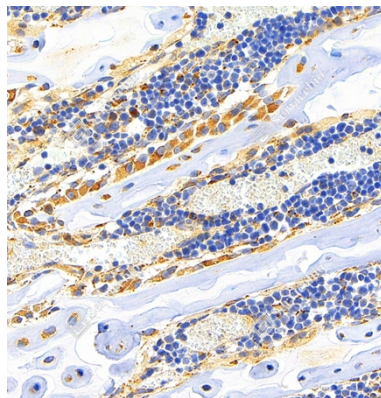

IHC analysis of RANKL(GB115705).
Sample: Rat bone, 4% PFA (G1101) 12-24h.
Antigen retrieval: Citrate buffer (pH 6.0) (G1201), pressure cooker 2 min.
Primary antibody: 1: 1000, 4°C overnight.
Secondary antibody: S-vision poly-HRP conjugated Goat Anti-Rabbit IgG (H+L), Ready to use(G1302), RT, 20min.

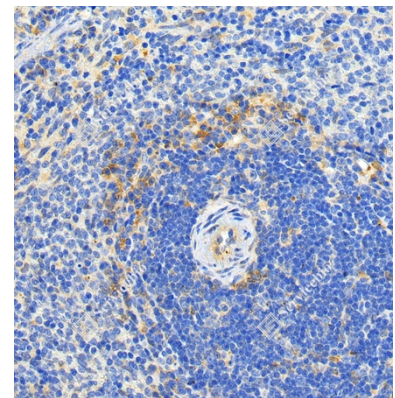

IHC analysis of RANKL(GB115705).
Sample: Rat spleen, 4% PFA (G1101) 12-24h.
Antigen retrieval: Citrate buffer (pH 6.0) (G1201), pressure cooker 2 min.
Primary antibody: 1: 1000, 4°C overnight.
Secondary antibody: S-vision poly-HRP conjugated Goat Anti-Rabbit IgG (H+L), Ready to use(G1302), RT, 20min.

For Research Use Only!

Ver. No.: V2.0-20260420