


Anti-TRAF1 Rabbit pAb



GB11595 -20°C

Product Information			
Protein full name	TNF receptor-associated factor 1		
Synonyms	TRAF1, EBI6, MGC, 10353, TNF receptor associated factor 1		
Immunogen	Recombinant protein corresponding to Mouse TRAF1		
Uniprot ID	P39428		
Isotype	IgG		
Purity	Affinity purification		
Predicted MW./Observed MW.	46 kDa / 48 kDa		
Applications	Species	Dilution	Positive Tissue
WB	Mouse, Rat	1: 500-1: 1000	liver, spleen, lung
Storage			
Storage	Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		
Storage Buffer	PBS with 0.15% ProClin300, 100 µg/mL BSA and 50% glycerol.		
background			
The protein encoded by this gene is a member of the TNF receptor (TNFR) associated factor (TRAF) protein family. TRAF proteins associate with, and mediate the signal transduction from various receptors of the TNFR superfamily. This protein and TRAF2 form a heterodimeric complex, which is required for TNF-alpha-mediated activation of MAPK8/JNK and NF-kappaB. The protein complex formed by this protein and TRAF2 also interacts with IAP, and thus mediates the anti-apoptotic signals from TNF receptors.			
Images			

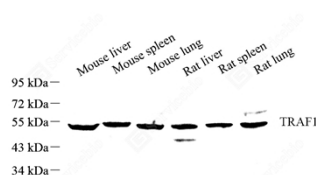

WB analysis of TRAF1(GB11595).
Sample: Protein treated by RIPA Lysis Buffer (G2002).
Blocking buffer: 3% Nonfat dry milk (GC310001) in TBST, RT, 1h.
Primary antibody: 1: 500, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 5000, RT, 1h.

For Research Use Only!

Ver. No.: V2.0-20260420