


Anti-TRAF3 Rabbit pAb




GB11722 -20°C

Product Information			
Protein full name	TNF receptor-associated factor 3		
Synonyms	TRAF3, Traf3, CAP-1, CAP1, CD40bp, CRAF1, IIAE5, LAP1, TNF receptor associated factor 3, RNF118		
Immunogen	Recombinant protein corresponding to Mouse TRAF3		
Uniprot ID	Q60803		
Isotype	IgG		
Purity	Affinity purification		
Subcellular location	Cytoplasm		
Applications	Species	Dilution	Positive Tissue
IHC/IF	Mouse	1: 500-1: 3000	hepatitis, spleen
Storage			
Storage	Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		
Storage Buffer	PBS with 0.15% ProClin300, 100 µg/mL BSA and 50% glycerol.		
background			
TNF receptor-associated factor is a member of the TNF receptor associated factor (TRAF) protein family. TRAF proteins associate with, and mediate the signal transduction from, members of the TNF receptor (TNFR) superfamily. This protein participates in the signal transduction of CD40, a TNFR family member important for the activation of the immune response. This protein is found to be a critical component of the lymphotoxin-beta receptor (LTbetaR) signaling complex, which induces NF-kappaB activation and cell death initiated by LTbeta ligation.			
Images			

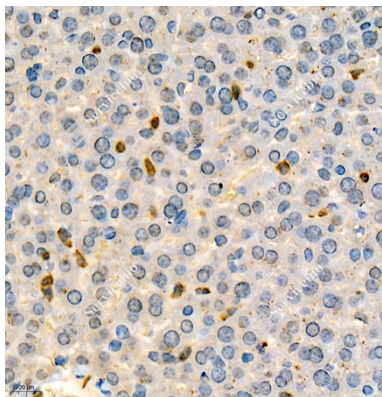

IHC analysis of TRAF3(GB11722).
Sample: Mouse hepatitis, 4% PFA (G1101) 12-24h.
Antigen retrieval: Citrate buffer (pH 6.0) (G1201), 98°C, 20 min.
Primary antibody: 1: 1500, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 200, RT, 1h.

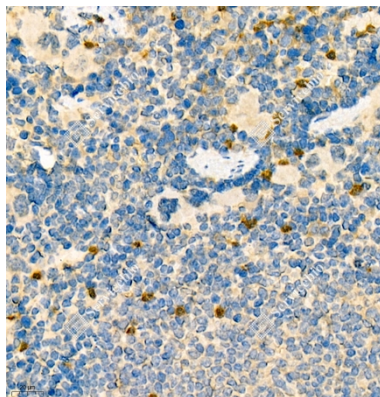

IHC analysis of TRAF3(GB11722).
Sample: Mouse spleen, 4% PFA (G1101) 12-24h.
Antigen retrieval: Citrate buffer (pH 6.0) (G1201), 98°C, 20 min.
Primary antibody: 1: 1500, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 200, RT, 1h.

For Research Use Only!

Ver. No.: V2.0-20260420