


Anti-MALT1 Rabbit pAb




GB11790 -20°C


Product Information			
Protein full name	Mucosa-associated lymphoid tissue lymphoma translocation protein 1 homolog		
Synonyms	Malt1, MALT1, IMD12, MLT, MLT1, PCASP1, MALT1 paracaspase		
Immunogen	Recombinant protein corresponding to Mouse MALT1		
Uniprot ID	Q2TBA3		
Isotype	IgG		
Purity	Affinity purification		
Subcellular location	Cytoplasm, Nucleus		
Applications	Species	Dilution	Positive Tissue
IHC/IF	Mouse, Rat	1: 250-1: 2000	lymph node, spleen, thymus
Storage			
Storage	Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		
Storage Buffer	PBS with 0.15% ProClin300, 100 µg/mL BSA and 50% glycerol.		
background			
Enhances BCL10-induced activation of NF-kappa-B. Involved in nuclear export of BCL10. Binds to TRAF6, inducing TRAF6 oligomerization and activation of its ligase activity. Has ubiquitin ligase activity. MALT1-dependent BCL10 cleavage plays an important role in T-cell antigen receptor-induced integrin adhesion.			
Images			

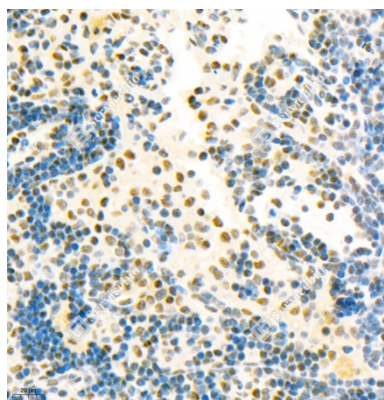

IHC analysis of MALT1(GB11790).
Sample: Mouse lymph node, 4% PFA (G1101) 12-24h.
Antigen retrieval: Citrate buffer (pH 6.0) (G1201), 98°C, 20 min.
Primary antibody: 1: 1000, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 200, RT, 1h.

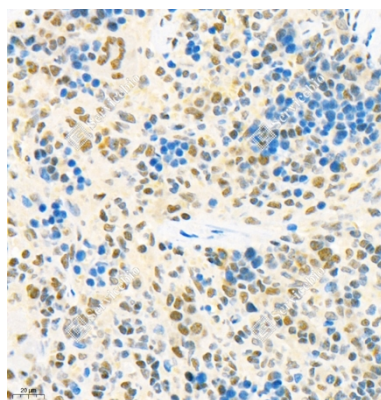

IHC analysis of MALT1(GB11790).
Sample: Rat lymph node, 4% PFA (G1101) 12-24h.
Antigen retrieval: Citrate buffer (pH 6.0) (G1201), 98°C, 20 min.
Primary antibody: 1: 1000, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 200, RT, 1h.

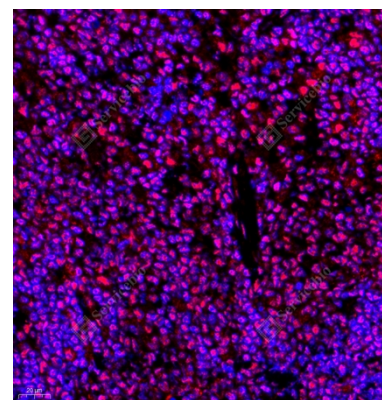

IF analysis of MALT1(GB11790)red.
Sample: Rat spleen, 4% PFA (G1101) 12-24h.
Antigen retrieval: Tris-EDTA buffer (pH 9.0) (G1203), 98°C, 20min.
Blocking buffer: 3% BSA (GC305010) in PBS, RT, 30min.
Primary antibody: 1: 500, 4°C overnight.
Secondary antibody: Cy3 conjugated Goat Anti-Rabbit IgG (H+L) (GB21303), 1: 300, RT, 1h.

For Research Use Only!

Ver. No.: V2.0-20260420