


Anti-FLIP Rabbit pAb



GB11817 -20°C

Product Information			
Protein full name	CASP8 and FADD-like apoptosis regulator		
Synonyms	CASH, Cflar, CFLAR, CASP8AP1, CLARP, Casper, FLAME, FLAME-1, FLAME1, FLIP, I-FLICE, MRIT, c-FLIP, c-FLIPL, c-FLIPR, c-FLIPS, CASP8 and FADD like apoptosis regulator		
Immunogen	Recombinant protein corresponding to Mouse FLIP		
Uniprot ID	O35732		
Isotype	IgG		
Purity	Affinity purification		
Predicted MW./Observed MW.	55 kDa / 55 kDa		
Applications	Species	Dilution	Positive Tissue
WB	Human, Mouse, Rat	1: 500-1: 1000	kidney, heart
Storage			
Storage	Store at -20°C for one year. Avoid repeated freeze/thaw cycles.		
Storage Buffer	PBS with 0.15% ProClin300, 100 µg/mL BSA and 50% glycerol.		
background			
CASP8 and FADD-like apoptosis regulator is a protein encoded by the CFLAR gene. It is a regulator of apoptosis and is structurally similar to caspase-8. However, the encoded protein lacks caspase activity and appears to be itself cleaved into two peptides by caspase-8. Several transcript variants encoding different isoforms have been found for this gene, and partial evidence for several more variants exists.			
Images			

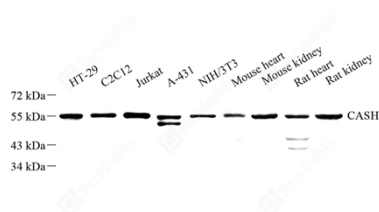

WB analysis of FLIP(GB11817).
Sample: Protein treated by RIPA Lysis Buffer (G2002).
Blocking buffer: 3% Nonfat dry milk (GC310001) in TBST, RT, 1h.
Primary antibody: 1: 500, 4°C overnight.
Secondary antibody: HRP conjugated Goat Anti-Rabbit IgG (H+L) (GB23303), 1: 5000, RT, 1h.

For Research Use Only!

Ver. No.: V2.0-20260420